

*t*rigon[®]


CALCANEAL Z WEDGE SYSTEM

The Trigon HA Stand-Alone Wedge Fixation System is a family of PEEK Optima HA Enhanced (HA PEEK) wedges with tantalum markers used for angular correction of small bones of the foot.

The wedges incorporate two screw-receiving holes, surface teeth, and an area to contain grafting material. The wedges are designed in rectangular, kidney, circular and oval shaped footprints in a range of sizes and in multiple thicknesses. The associated 2.5mm diameter titanium screws are designed in lengths of 10 to 30mm.

When used with the provided screw fixation the Trigon wedges may be used with or without ancillary plating, except for the Lapidus Wedges which will require ancillary fixation not provided in the Trigon System. When used without the provided screws, Trigon wedges are intended for use with ancillary fixation.

- Fusion Grade PEEK is a unique compound of PEEK (80%) and hydroxyapatite (20%) for radiolucent stable fixation with osteoconductive and osteoinductive properties.
- A cut guide with depth markings is included to trim nails as needed.
- Tantalum markers at the distal end of each nail provide a depth indicator on imaging.
- Cannulated for k-wire guidance for enhanced surgical accuracy.
- Nails are impacted across application sites to provide solid fixation.

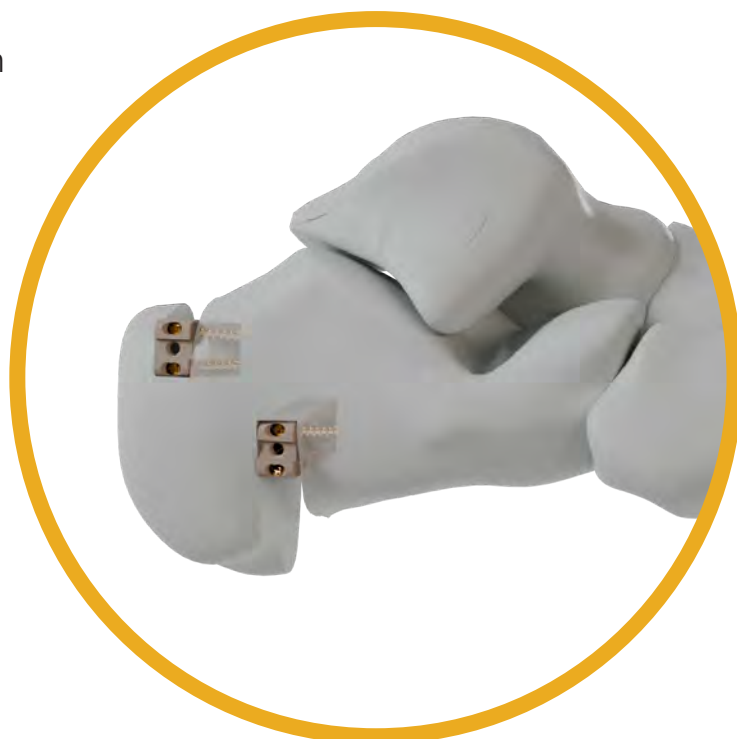

Ordering Information



Wedges

Part Number	Description
CZP-1020-05	PEEK Calcaneal Z Evans Wedge 5mm
CZP-1020-75	PEEK Calcaneal Z Evans Wedge 7.5mm
CZP-1020-10	PEEK Calcaneal Z Evans Wedge 10mm
CZP-1535-05	PEEK Calcaneal Z Small Tuber Wedge 5mm
CZP-1535-75	PEEK Calcaneal Z Small Tuber Wedge 7.5mm
CZP-1535-10	PEEK Calcaneal Z Small Tuber Wedge 10mm
CZP-2025-05	PEEK Calcaneal Z Large Tuber Wedge 5mm
CZP-2025-75	PEEK Calcaneal Z Large Tuber Wedge 7.5mm
CZP-2025-10	PEEK Calcaneal Z Large Tuber Wedge 10mm

Trials

Part Number	Description
CZP-1020T-05	PEEK Calcaneal Z Evans Trial Wedge 5mm
CZP-1020T-75	PEEK Calcaneal Z Evans Trial Wedge 7.5mm
CZP-1020T-10	PEEK Calcaneal Z Evans Trial Wedge 10mm
CZP-1535T-05	PEEK Calcaneal Z Small Tuber Trial Wedge 5mm
CZP-1535T-75	PEEK Calcaneal Z Small Tuber Trial Wedge 7.5mm
CZP-1535T-10	PEEK Calcaneal Z Small Tuber Trial Wedge 10mm
CZP-2025T-05	PEEK Calcaneal Z Large Tuber Trial Wedge 5mm
CZP-2025T-75	PEEK Calcaneal Z Large Tuber Trial Wedge 7.5mm
CZP-2025T-10	PEEK Calcaneal Z Large Tuber Trial Wedge 10mm

Instruments

Part Number	Description
WD-1000T-ZSHFT	Calcaneal Z Wedge Inserter Shaft
WD-1000T-SREW	Trigon Screw Inserter
WD-1000T-TRI	Trigon Trial Inserter
WD-1000T-DRL	Trigon 2.0mm Drill
CSRW-1000T-352	Healix Compression Screw K-Wires (2.4mm x 150mm)
HAN-1000T-AOSM	Handle
WD-1000T-CZER	Calcaneal Z Evans Cut Guide Right
WD-1000T-CZEL	Calcaneal Z Evans Cut Guide Left
WD-1000T-CZTR	Calcaneal Z Tuber Cut Guide Right
WD-1000T-CZTL	Calcaneal Z Tuber Cut Guide Left
WD-1001T-ZWSL	Calcaneal Z Wedge Crossed Screw Sleeve
WD-1000T-CZGH	Calcaneal Z Cut Guide Handle
WD-1001T-SDWSL	Calcaneal Z Wedge Convergent Screw Sleeve
	Hintermann
112-01	Bone Graft Harvester (6mm)
112-02	Bone Graft Harvester (8mm)
112-04	Bone Graft Harvester (10mm)




Nvision Biomedical Technologies

12449 Silicon Drive
San Antonio, TX 78249

☎ 210-545-3713 ext. 112

✉ orders@nvisionbiomed.com

🌐 www.nvisionbiomed.com

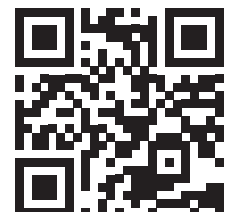

Scan to visit
our website.