

*t*rigon[®]


LESSER METATARSAL WEDGE SYSTEM

The Trigon Lesser Metatarsal Wedge System is made of Fusion Grade PEEK, where hydroxyapatite (HA) is fully integrated throughout the PEEK.

The Trigon Lesser Metatarsal Wedge System is indicated for a metatarsal-cuneiform lengthening arthrodesis of the second and/or third TMT joint(s). Both single and double joint options are provided, in parallel or angular correction. The wedge is available in a variety of sizes that accommodate a large patient population, allowing precise and reliable correction.

The Trigon Lesser Metatarsal Wedge System can be combined with ancillary fixation via the staples provided. A targeting guide is provided for screw fixation if preferred.

- Fusion Grade PEEK is a unique compound of PEEK (80%) and hydroxyapatite (20%) for radiolucent stable fixation with osteoconductive and osteoinductive properties.
- A cut guide with depth markings is included to trim nails as needed.
- Tantalum markers at the distal end of each nail provide a depth indicator on imaging.
- Cannulated for k-wire guidance for enhanced surgical accuracy.
- Nails are impacted across application sites to provide solid fixation.

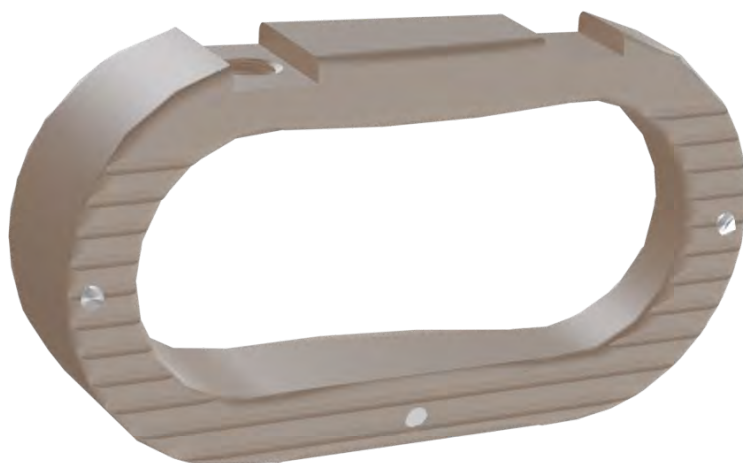

Ordering Information

Wedges SINGLE JOINT

Part Number	Description
LMP-1414-0600	6mm Wedge Parallel PEEK Small Single Joint Lesser Metatarsal
LMP-1414-0806	8mm Wedge 6° Transverse PEEK Small Single Joint Lesser Metatarsal
LMP-1414-0809	8mm Wedge 9° Transverse PEEK Small Single Joint Lesser Metatarsal
LMP-1715-0806	8mm Wedge 6° Transverse PEEK Large Single Joint Lesser Metatarsal
LMP-1715-0809	8mm Wedge 9° Transverse PEEK Large Single Joint Lesser Metatarsal
LMP-1428-0600	6mm Wedge Parallel PEEK Small Multi Joint Lesser Metatarsal
LMP-1428-0806	8mm Wedge 6° Transverse PEEK Small Multi Joint Lesser Metatarsal
LMP-1428-0809	8mm Wedge 9° Transverse PEEK Small Multi Joint Lesser Metatarsal
LMP-1731-0806	8mm Wedge 6° Transverse PEEK Large Multi Joint Lesser Metatarsal
LMP-1731-0809	8mm Wedge 9° Transverse PEEK Large Multi Joint Lesser Metatarsal

Staples

Part Number	Description
STK-20-14-14	20mm Bridge Staple, 14mm x 14mm Legs
STK-20-18-18	20mm Bridge Staple, 18mm x 18mm Legs
STK-20-20-20	20mm Bridge Staple, 20mm x 20mm Legs

Trials

Part Number	Description
LM-1414T-0600	6mm Trial Parallel PEEK Small Single Joint Lesser Metatarsal
LM-1414T-0806	8mm Trial 6° Transverse PEEK Small Single Joint Lesser Metatarsal
LM-1414T-0809	8mm Trial 9° Transverse PEEK Small Single Joint Lesser Metatarsal
LM-1715T-0806	8mm Trial 6° Transverse PEEK Large Single Joint Lesser Metatarsal
LM-1715T-0809	8mm Trial 9° Transverse PEEK Large Single Joint Lesser Metatarsal
LM-1428T-0600	6mm Trial Parallel PEEK Small Multi Joint Lesser Metatarsal
LM-1428T-0806	8mm Trial 6° Transverse PEEK Small Multi Joint Lesser Metatarsal
LM-1428T-0809	8mm Trial 9° Transverse PEEK Small Multi Joint Lesser Metatarsal
LM-1731T-0806	8mm Trial 6° Transverse PEEK Large Multi Joint Lesser Metatarsal
LM-1731T-0809	8mm Trial 9° Transverse PEEK Large Multi Joint Lesser Metatarsal

Instruments

Part Number	Description
HAN-1000T-AOSM	Handle
WD-1000T-LMCG06	Lesser Metatarsal 6mm Cut Guide
WD-1000T-LMCG08	Lesser Metatarsal 8mm Cut Guide
WD-1000T-LMSDG	Lesser Metatarsal Wedge Staple Drill Guide
WD-1000T-MCGH	Lesser Metatarsal/Midfoot Cut Guide Handle
WD-1000T-TRI	Trigon Trial Insertor
STK-1000T-011	2.5mm Drill
STK-1000T-DIS1	Staple Distractor
STK-1000T-TMP	Tamp

