

*radian*TM

MIS BUNION CORRECTION SYSTEM

The RadianTM MIS Bunion System is indicated for fixation of osteotomies and correction of hallux valgus. The construct consists of a titanium alloy plate and four locking screws that provide fixation for the correction of a first metatarsal bunion. The design allows for medial insertion of the two distal screws and intramedullary insertion of the two proximal screws.

The RadianTM MIS Bunion System allows the plate to be attached to the metatarsal head prior to making the osteotomy. This allows for control of the capital fragment at all stages of the procedure.

The plate is provided in two configurations, a gold plate for IM shifts greater than 10 deg, a blue plate for IM shift less than 10 deg.



Ordering Information

Plates

Part Number	Description
BUN-PLT-05	5.0mm Correction L-Bunion Plate
BUN-PLT-03	5.0mm Correction L-Bunion Plate 30 deg

Screws

Part Number	Description
TB-30-14L	3.0mmD, 14mmL Headed, Locking Screw
TB-30-16L	3.0mmD, 16mmL Headed, Locking Screw
TB-30-18L	3.0mmD, 18mmL Headed, Locking Screw
TB-30-20L	3.0mmD, 20mmL Headed, Locking Screw
TB-30-22L	3.0mmD, 22mmL Headed, Locking Screw
TB-30-24L	3.0mmD, 24mmL Headed, Locking Screw
TB-30-26L	3.0mmD, 26mmL Headed, Locking Screw
TB-30-28L	3.0mmD, 28mmL Headed, Locking Screw
TB-30-30L	3.0mmD, 30mmL Headed, Locking Screw

Instruments


Part Number	Description
BUN-1000T-DGID	Drill Guide
BUN-1000T-DR30	3.0mm Diameter Drill
BUN-1000T-OPIN	Olive Pin
TB-1000T-ISRT	Screw Inserter
BUN-1001T-DPG	Depth Gauge
HAN-1000T-AOSM	Handle
BUN-1000T-SPIN	Small Olive Pin
CSRW-1000T-352	2.4mm x 150mm K-Wire
BUN-1000T-SGID	Drill and Screw Guide
BUN-1000T-SGHS	Drill and Screw Guide Holding Screw
BUN-1000T-KDG	K-wire Depth Gauge
BUN-1000T-SG30	Drill and Screw Guide 30°
BUN-1000T-SGHS	Drill and Screw Guide Holding Screw



NVISION Biomedical Technologies

12449 Silicon Drive
San Antonio, TX 78249

210-545-3713 ext. 112
orders@nvisionbiomed.com
www.nvisionbiomed.com

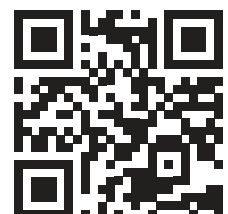

Scan to visit
our website.